

Instructions for Testing Residue for IV Drug Users

If you inject heroin or other drugs, you should test your drugs every time you inject.

The easiest method is to test the residue from your spoon or cooker:



1. After preparing your shot, set the needle aside and wait to inject.


2. Add about 1 ml (1/4 of a teaspoon) of clean water into the spoon or cooker.

3. Swish the water around inside the spoon or cooker.

4. Hold the blue end of the test strip and insert the other end into the liquid (no higher than the blue line) for about 15 seconds.

5. Set the strip down on a flat surface and wait one to two minutes.

6. Read the results.

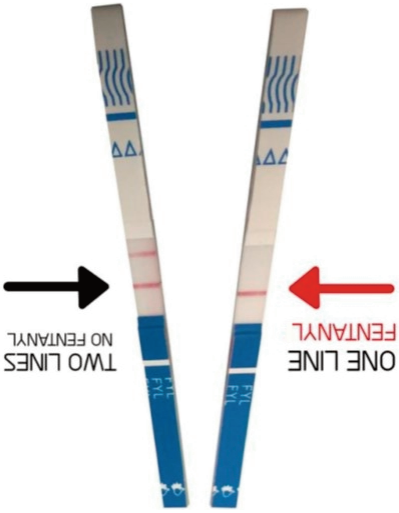

INSTRUCTIONS ON HOW TO USE FENTANYL TEST STRIPS



Important Information About Fentanyl Test Strips

Fentanyl test strips are for harm reduction purposes only. They cannot tell you how much fentanyl is present, and they cannot detect every fentanyl analog or other synthetic opiates. A negative test does not mean the drugs being tested are 100% safe. All street drug use contains inherent risks.

How to Correctly Read the Test Results

One Red Line = Fentanyl

Two Red Lines = No Fentanyl

No Lines or One Red Line at the Bottom =

The test didn't work!

Safety Tips

- Don't use alone.
- Always carry Narcan®.
- Call **911** in case of an overdose.



Instructions for Testing Cocaine and Other Non-Injected Drugs (Best Method)

1. Place all the drugs you intend to use into a small glass or ceramic cup. (Crush tablets and

crystals into a powder.)

2. Add water.
 - Add 2 teaspoons of water per 100 mg of powder (~10 mg/ml).

3. Stir thoroughly with a clean spoon.

4. Set the spoon down.

5. Use another spoon to fill the residue spoon with clean water.

6. Hold the blue end of the test strip and insert the other end into the residue liquid in the spoon for 15 seconds.

5. Set the strip down on a flat surface and wait one to two minutes.

6. Read the results.

CAUTION: If you are testing meth or MDMA, add 1 teaspoon of water (~5 ml) for each 10 mg of crystal or powder to avoid a false positive.

Instructions for Testing the Residue Inside a Baggie

CAUTION: This is not the best method, because it is best to test all that you intend to consume each time you use (as described in the "Best Method" to the right).

NOTE: If you are testing cocaine that is pressed into a rock, grind it or crush it up and put it in a baggie before following the steps below.

1. Shake your baggie, then empty the powder or crystals inside your baggie into a second baggie.

2. Fill the first baggie with 1/2 teaspoon of clean water and swish it around to dissolve the residue.

- If testing meth or MDMA, you may need to add more water depending on how much residue is stuck to the inside of the baggie.

3. Pour the liquid-residue into a clean glass and add half a cup of water to dilute the residue.

4. Hold the blue end of the test strip and insert the other end into the residue liquid in the glass for 15 seconds.

5. Set the strip down on a flat surface and wait one to two minutes.

6. Read the results.

Please help us by taking a private survey about fentanyl test strips:

Visit this web page

www.surveymonkey.com/r/5FN59MZ

Or use this QR code

